

Technical Specifications

Wavelength range	190 – 1,100 nm
Wavelength scan range	200 – 950 nm
Measure time for full scan range	Less than 5 seconds
Wavelength reproducibility	< ± 0.2 nm
Wavelength accuracy	± 2 nm
Bandwidth	Better than 5 nm
Stray light	< 0.5% at 220 nm and 340 nm using NaNO ₂
Photometric range	-0.3 – 2.499 A 0-199% T
Photometric reproducibility	± 0.003 A (0 to 0.5 A) ± 0.007 A (0.5-1.0 A) @ 260 nm
Photometric accuracy	± 0.005 A or $\pm 1\%$ of the reading, whichever is the greater
Zero stability	± 0.003 A/hour after 20 min warm up @ 340 nm
Noise	0.002 A rms at 0 A @ 260 nm 0.005 A (pk to pk) at 0 A @ 260 nm
Optical arrangement	Dual channel Czerny Turner with flat grating, 1024 pixel CCD array, concave mirrors.
Lamp	Xenon flash lamp
Lifetime	10 ⁹ flashes, up to 10 years
Warranty	1 year
Performance verification	Auto diagnostics when switched on
Cell types	15 mm centre height, outside dimension 12.5 mm x 12.5 mm
Cuvette storage	capacity for eight 10 mm cells
Photometric mode	Abs, %T, concentration, scan, ratio, multi wavelength, kinetics in Δ Abs x factor/min
Method storage	Up to 90 methods in user methods
Built-in methods	Nucleid acid, microarray (labeling efficiency), protein and cell density
Display formats	320 x 240 pixels
Size	140 mm x 275 mm x 380 mm
Weight	< 4.5 kg
Operating voltage	90-250 V, 50/60 Hz, Max 30 VA
Input / Output ports	USB or Bluetooth for connection to a PC for direct data download for spreadsheet calculations, printout and data storage

Features and specifications are subject to change without notice.



Worldwide
Headquarter

IMPLEN GmbH
Wehrlestr. 33
81679 München
Germany

Tel.: +49 89 99100583
Fax: +49 89 21758349
Email: info@implen.de
www.implen.de

UK Office
Implen UK Ltd.
Cumberland House
24-28 Baxter Avenue
Southend on Sea
Essex, SS2 6HZ
Great Britain

Tel.: +44 1702 335588
Fax: +44 1702 334262
Email: info@implen.co.uk
www.implen.co.uk

The NanoPhotometer™



The NanoPhotometer™

One spectrophotometer for all your needs

All kinds of UV/Vis photometric applications in molecular biology, biochemistry and microbiology featuring

- ultra low sample volume: starting from 0.7 μ l
- measurements without cuvettes plus standard cuvette option

- predefined methods for protein and nucleic acid measurements
- full spectrum scan and kinetic methods
- full developed software on-board and direct data download to PC
- high reproducibility

SMALL & FAST

Print data via built-in printer, store data via USB or Bluetooth on your PC

READY TO USE

Low bench space requirement, computer independent. Built-in computer with backlit LCD display for fast and easy data access

ALL IN ONE



Multi purpose sample holder for sub-microliter, standard cuvette and fiber optic applications

ULTRA LOW VOLUME, NO SAMPLE DILUTION

LabelGuard™ microliter cell device for sub-microliter measurements included within the system

EASY TO CONTROL

Rugged key pad for fast method selection and easy data entry

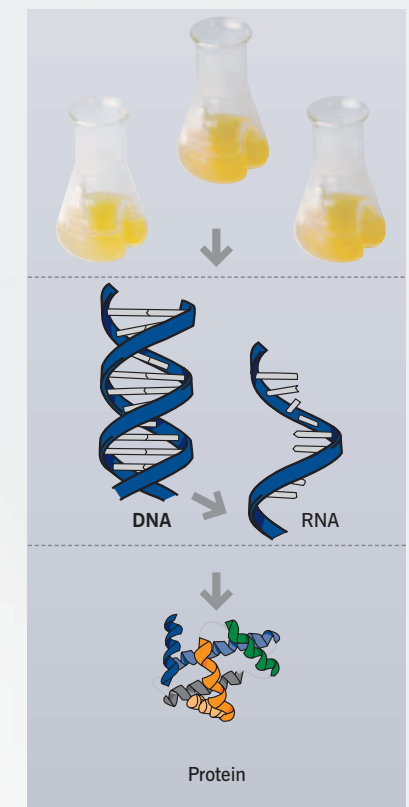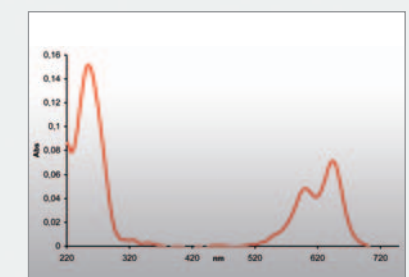

VOLUME

WAVELENGTH

CELL DENSITY

NUCLEIC ACIDS

PROTEIN (BSA)

RANGE

3,500 μ l

0.7 μ l

1,100 nm

190 nm

8-10⁸ cells/ml

1 cell/ml

4,250 ng/μl

15 ng/μl

100 mg/ml

0.2 mg/ml